

For Laboratory Use Only

PRODUCT NAME	Genistein	Structure
PART NUMBER	00007082	
STANDARD TYPE	Primary (P)	HO
LOT NUMBER	00007082-550	
REPORT NUMBER	CDXA-RSS-8897-00	
SAMPLE NUMBER	CDXP-20-00193	ОН
DATE OF SAMPLE	03/30/2020	∽он
DATE OF REPORT	04/30/2020	
CHEMICAL NAMES	4',5,7-Trihydroxyisoflavone; Genisteol	; Genisterin; Prunetol; Sophoricol

	/
CHEMICAL FORMULA	$C_{15}H_{10}O_5$
MOLECULAR WEIGHT (MW)	270.24
CHEMICAL FAMILY	Flavonoids
CAS NUMBER	[446-72-0]
EC#(EINECS)	207-174-9
RTECS	NR2392000

ANALYTICAL RESULTS

TEST	METHOD	SPECIFICATION	RESULT
Adjusted Purity	NA	NA	98.9%
HPLC	0.700.10.2.METH22	NA	99.3%
NMR	0.700.12.4	Conforms to structure	Conforms to structure
Mass Spectrum	0.700.12.27	Conforms	Conforms
Residual Solvent*	USP <467>	NA	Methanol – 0.4%
Water	0.700.12.37	NA	ND
Appearance	NA	NA	White powder

*Testing performed at an ISO 17025 accredited subcontracted laboratory. ADJUSTED PURITY: 98.9% IS BASED ON (100 - 0.4 SOLVENTS - 0.0 WATER) X 99.3% HPLC

STORAGE CONDITIONS

STORAGE -20 °C in a dry place.

EXPIRATION DATE 03/2025 under the above conditions.

This document is the property of ChromaDex, Inc. and/or its relevant affiliates and contains confidential and proprietary material for the sole use of the intended recipient(s). Any review, use, distribution or disclosure by or to others is strictly prohibited.

10005 MUIRLANDS BLVD., SUITE G, IRVINE, CA 92618 | T: +1 949-419-0288 | F: +1 949-419-0294 | WWW.CHROMADEX.COM LABORATORY ADDRESS: 1751 S. FORDHAM ST., SUITE 350, LONGMONT, CO 80503

India Contact:


306, Aggarwal City Mall, Opposite M2K Pitampura, Delhi – 110034 (INDIA). Ph: +91-11-42208000, 42208111, 42208222, Mobile: +91-9810521400, Fax: +91-11-42208444 Email: customerservice@lifetechindia.com



ANALYTICAL CONDITIONS

INSTRUMENT	AGILENT 1260 HPLC UV-VIS (DAD); AB-SCIEX 6600 TRIPLE TOF
COLUMN	Phenomenex Luna C18(2) 150 x 2.00 mm, 5 μ m particle size
MOBILE PHASE	A - 0.1% Trifluoroacetic acid in Ultrapure water, B - Acetonitrile; Isocratic 5% B for 5 minutes, then increasing to 95% B over 20 minutes. Hold at 95% B for 15 minutes, then decreasing to 5% B over 1 minute.
COLUMN TEMP.	0° C
FLOW RATE	0.4 mL/minute
INJECTION VOL.	0.2 μL
INJECTION CONC.	1.1 mg/mL in Methanol
DETECTION	205 ± 4 nm

HPLC CHROMATOGRAM OF GENISTEIN (CDXP-20-00193)

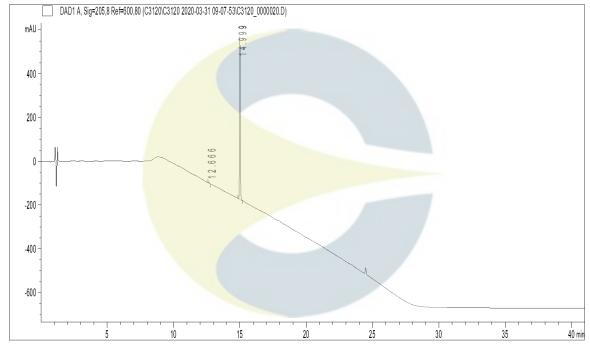


This document is the property of ChromaDex, Inc. and/or its relevant affiliates and contains confidential and proprietary material for the sole use of the intended recipient(s). Any review, use, distribution or disclosure by or to others is strictly prohibited.

India Contact: Life Technologies (India) Pvt. Ltd. 306, Aggarwal City Mall, Opposite M2K Pitampura, Delhi – 110034 (INDIA). Ph: +91-11-42208000, 42208111, 42208222, Mobile: +91-9810521400, Fax: +91-11-42208444 Email: customerservice@lifetechindia.com Website: www.lifetechindia.com




UV SPECTRUM AT RT 15.0 MINUTES

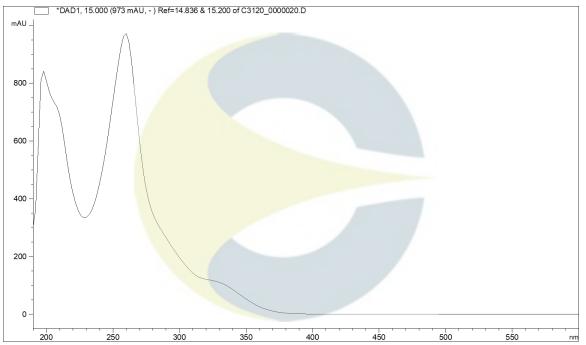

HPLC CHROMATOGRAM OF DILUENT BLANK

India Contact:

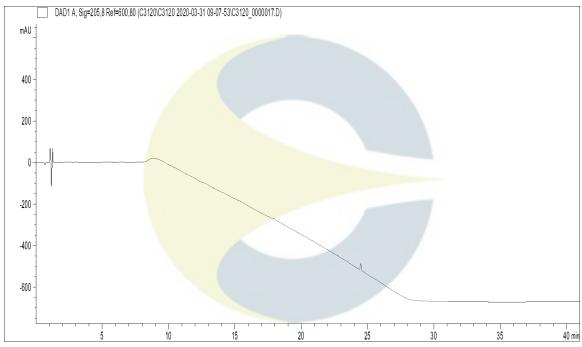

This document is the property of ChromaDex, Inc. and/or its relevant affiliates and contains confidential and proprietary material for the sole use of the intended recipient(s). Any review, use, distribution or disclosure by or to others is strictly prohibited.

10005 MUIRLANDS BLVD., SUITE G, IRVINE, CA 92618 | T: +1 <u>949-419-0288</u> | F: +1 <u>949-419-0294</u> | WWW.CHROMADEX.COM LABORATORY ADDRESS: 1751 S. FORDHAM ST., SUITE 350, LONGMONT, CO 80503

Life Technologies (India) Pvt. Ltd. 306, Aggarwal City Mall, Opposite M2K Pitampura, Delhi – 110034 (INDIA). Ph: +91-11-42208000, 42208111, 42208222, Mobile: +91-9810521400, Fax: +91-11-42208444 Email: customerservice@lifetechindia.com Website: www.lifetechindia.com


For Laboratory Use Only




India Contact:

For Laboratory Use Only

MASS SPECTRUM OF GENISTEIN (CDXP-20-00193)

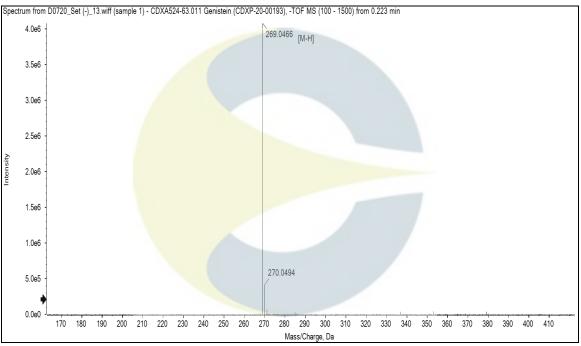

¹H-NMR SPECTRUM OF GENISTEIN (CDXP-20-00193) IN DMSO-d₆

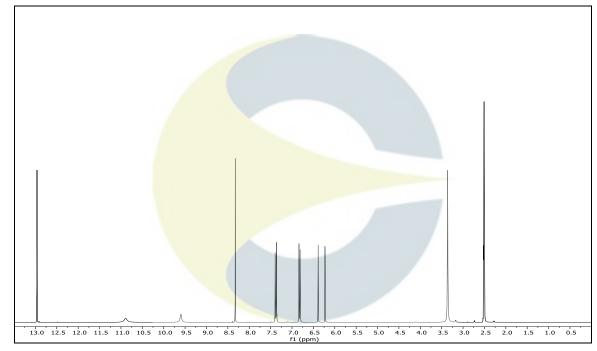

This document is the property of ChromaDex, Inc. and/or its relevant affiliates and contains confidential and proprietary material for the sole use of the intended recipient(s). Any review, use, distribution or disclosure by or to others is strictly prohibited.

10005 MUIRLANDS BLVD., SUITE G, IRVINE, CA 92618 | T: +1 <u>949-419-0288</u> | F: +1 <u>949-419-0294</u> | WWW.CHROMADEX.COM LABORATORY ADDRESS: 1751 S. FORDHAM ST., SUITE 350, LONGMONT, CO 80503


Life Technologies (India) Pvt. Ltd. 306, Aggarwal City Mall, Opposite M2K Pitampura, Delhi – 110034 (INDIA). Ph: +91-11-42208000, 42208111, 42208222, Mobile: +91-9810521400, Fax: +91-11-42208444 Email: customerservice@lifetechindia.com Website: www.lifetechindia.com



India Contact:

For Laboratory Use Only

¹³C-NMR SPECTRUM GENISTEIN (CDXP-20-00193) IN DMSO-d₆

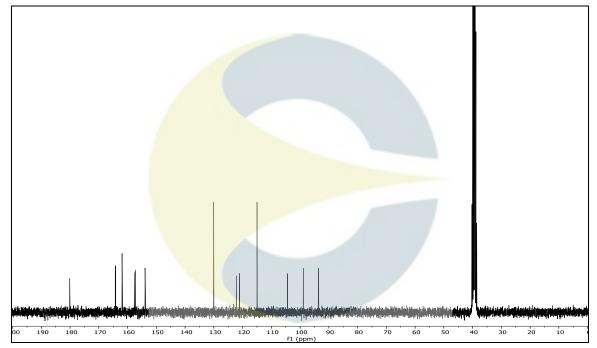


This document is the property of ChromaDex, Inc. and/or its relevant affiliates and contains confidential and proprietary material for the sole use of the intended recipient(s). Any review, use, distribution or disclosure by or to others is strictly prohibited.

Life Technologies (India) Pvt. Ltd. 306, Aggarwal City Mall, Opposite M2K Pitampura, Delhi – 110034 (INDIA). Ph: +91-11-42208000, 42208111, 42208222, Mobile: +91-9810521400, Fax: +91-11-42208444 Email: customerservice@lifetechindia.com Website: www.lifetechindia.com



For Laboratory Use Only

COSY-NMR SPECTRUM OF GENISTEIN (CDXP-20-00193) IN DMSO-d6

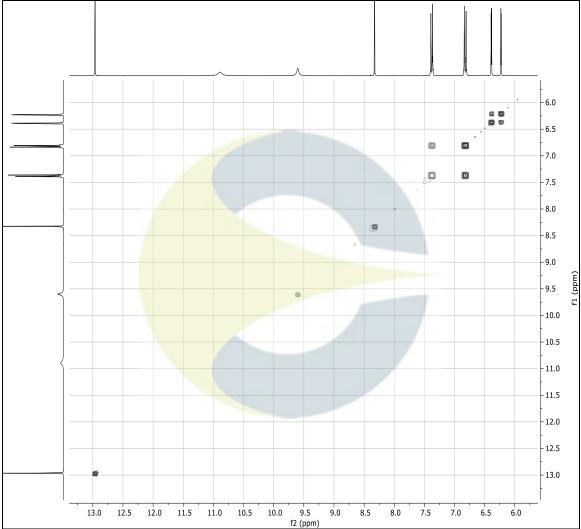

REVISION HISTORY

India Contact:

<u>Revision History</u>	Date of Revision	Document/Changes
00	04/30/2020	New report

This document is the property of ChromaDex, Inc. and/or its relevant affiliates and contains confidential and proprietary material for the sole use of the intended recipient(s). Any review, use, distribution or disclosure by or to others is strictly prohibited.

10005 MUIRLANDS BLVD., SUITE G, IRVINE, CA 92618 | T: +1 <u>949-419-0288</u> | F: +1 <u>949-419-0294</u> | WWW.CHROMADEX.COM LABORATORY ADDRESS: 1751 S. FORDHAM ST., SUITE 350, LONGMONT, CO 80503

Life Technologies (India) Pvt. Ltd. 306, Aggarwal City Mall, Opposite M2K Pitampura, Delhi – 110034 (INDIA). Ph: +91-11-42208000, 42208111, 42208222, Mobile: +91-9810521400, Fax: +91-11-42208444 Email: customerservice@lifetechindia.com Website: www.lifetechindia.com