


Dihydrorhodamine 123

Cat. No.:	HY-101894
CAS No.:	109244-58-8
Molecular Formula:	C ₂₁ H ₁₈ N ₂ O ₃
Molecular Weight:	346.38
Target:	Others
Pathway:	Others
Storage:	-20°C, protect from light

* The compound is unstable in solutions, freshly prepared is recommended.

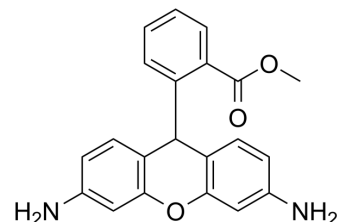

SOLVENT & SOLUBILITY

In Vitro

DMSO : 100 mg/mL (288.70 mM; Need ultrasonic)

Concentration	Mass			
	1 mg	5 mg	10 mg	
1 mM	2.8870 mL	14.4350 mL	28.8700 mL	
5 mM	0.5774 mL	2.8870 mL	5.7740 mL	
10 mM	0.2887 mL	1.4435 mL	2.8870 mL	

Please refer to the solubility information to select the appropriate solvent.

BIOLOGICAL ACTIVITY

Description

Dihydrorhodamine 123 (DHR 123) is a fluorescent probe (λ_{ex} =488 nm, λ_{em} =525 nm)^[1].

In Vitro

In the presence of 10 μ M Dihydrorhodamine 123 (DHR 123) the stimulation of the neutrophil NADPH oxidase by the addition of 50 nM phorbol 12-myristate 13-acetat (PMA) results in an increase in the rate of rhodamine generation. The fluorescent intensity of the cells, in the presence of 10 μ M Dihydrorhodamine 123, increases with time following the addition of 50 nM PMA. In the presence of 10 μ M Dihydrorhodamine 123, induced HL60 cells show a sustained increase in fluorescence following the addition of 50 nM PMA^[1].

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

PROTOCOL

Cell Assay ^[1]

The HL60 cells are incubated at 6×10^6 cells/mL in Krebs-Ringer buffer at 37°C containing 10 μ M Dihydrorhodamine 123 (DHR). The generation of O₂⁻ is initiated by the addition of 50 nM phorbol 12-myristate 13-acetat (PMA) and the progress of the generation of rhodamine 123 is monitored in 50- μ L aliquots (3×10^5 cells) diluted tenfold before analysis. The uninduced HL60 cells are loaded with 5 μ M carboxy SNARF-1 AM acetate (SNARF-AM) in the Na⁺ medium for 10 min at 37°C and washed

by centrifugation and resuspension to remove unhydrolysed SNARF ester^[1].
MCE has not independently confirmed the accuracy of these methods. They are for reference only.

CUSTOMER VALIDATION

- Part Fibre Toxicol. 2022 Mar 29;19(1):24.
- Oxid Med Cell Longev. 10 Jun 2022.
- J Mol Cell Cardiol. 2021 Jul 2;S0022-2828(21)00135-8.

See more customer validations on www.MedChemExpress.com

REFERENCES

[1]. Lydia M. Henderson et al. Dihydrorhodamine 123: a fluorescent probe for superoxide generation? Eur.J.Biochem. 217, 973-980.

Caution: Product has not been fully validated for medical applications. For research use only.

India Contact:
Life Technologies (India) Pvt. Ltd.
306, Aggarwal City Mall, Opposite M2K Pitampura, Delhi – 110034 (INDIA). Ph: +91-11-42208000, 42208111, 42208222, Mobile: +91-9810521400, Fax: +91-11-42208444
Email: customerservice@lifetechindia.com Website: www.lifetechindia.com