

## Canine N-terminal pro-brain natriuretic peptide (NT-proBNP) ELISA Kit

to determine Canine N-terminal pro-brain natriuretic peptide in Serum, Blood Plasma, Saliva, Urine, And Other Related Tissue Liquid Samples.

### INSTRUCTION MANUAL

FOR ELISA KIT No: LTaC9010EA



Arsh Biotech Pvt. Ltd.

308, Aggarwal City Mall, Road No.44,  
Pitampura, Delhi-110034, India

Mobile: +91-98105-21400 | Fax: +91-11-42208444

[info@arshbiotech.com](mailto:info@arshbiotech.com)

### INTENDED USE

This kit is used to assay the Canine N-terminal pro-brain natriuretic peptide (NT-proBNP) in the sample of serum, blood plasma, saliva, urine, and other related biological liquid. For *in vitro* use only.

MANUAL VERSION 1.02

### BACKGROUND

The N-terminal prohormone of brain natriuretic peptide (NT-proBNP or BNP) is a prohormone with a 76 amino acid N-terminal inactive protein that is cleaved from the molecule to release brain natriuretic peptide. Both BNP and NT-proBNP levels in the blood are used for screening, diagnosis of acute congestive heart failure (CHF) and may be useful to establish prognosis in heart failure, as both markers are typically higher in patients with worse outcome. The plasma concentrations of both BNP and NT-proBNP are also typically increased in patients with asymptomatic or symptomatic left ventricular dysfunction and is associated with coronary artery disease and myocardial ischemia.

### PRINCIPLE

Canine N-terminal pro-brain natriuretic peptide (NT-proBNP) ELISA Kit employs a two-site sandwich ELISA to quantitate NT-proBNP in samples. An antibody specific for NT-proBNP has been pre-coated onto a microplate. Standards and samples are pipetted into the wells and any NT-proBNP present is bound by the immobilized antibody. After removing any unbound substances, HRP-Conjugate Canine NT-proBNP detection antibody is added to the wells. Following a wash to remove any unbound HRP reagent, a Chromogen solution is added to the wells and color develops in proportion to the amount of NT-proBNP bound in the initial step. The color development is stopped and the intensity of the color is measured.

### ASSAY RESTRICTIONS

- Do not mix or substitute reagents with those from other lots or sources.
- It is important that the calibrator diluent selected for the standard curve be consistent with the samples being assayed.

- If samples generate values higher than the highest standard, dilute the samples with the appropriate calibrator diluent and repeat the assay.
- Any variation in standard diluent, operator, pipetting technique, washing technique, incubation time or temperature, and kit age can cause variation in binding.
- This assay is designed to eliminate interference by other factors present in biological samples. Until all factors have been tested in the ELISA Immunoassay, the possibility of interference cannot be excluded.

### MATERIALS SUPPLIED IN THIS KIT

**Canine NT-proBNP microplate:** 96 well polystyrene microplates (8 strips of 12 wells) coated with the antibody specific for Canine NT-proBNP.

**Canine NT-proBNP standard:** Canine NT-proBNP in a buffered protein base with preservatives, liquid.

**Standard diluent:** Diluent solution for reconstituted standard.

**Sample diluent:** Diluent solution for reconstituted samples.

**HRP-Conjugate Canine NT-proBNP detection antibody:** Antibody specific for Canine NT-proBNP, liquid.

**Chromogen solution A:** liquid.

**Chromogen solution B:** liquid.

**Stop solution:** liquid.

**Wash buffer:** 30X liquid.

Plate covers.

### MATERIALS REQUIRED BUT NOT SUPPLIED

1. 37 °C incubator
2. Microplate reader
3. Precision pipettes and tips
4. Distilled water
5. Disposable tubes for sample dilution
6. Absorbent paper

### STORAGE CONDITIONS

Store kit reagents at 2-8 °C for 12 months. Immediately after use remaining reagents should be returned to cold storage at 4°C. Recommended storage instruction for opened/reconstituted kit components are listed below.

Kit components	Storage conditions for opened components
Canine NT-proBNP microplate	Return unused wells to the foil pouch containing the desiccant pack. Reseal along entire edge of the zip-seal. May be stored for up to 1 month at 2-8 °C post opening.
Canine NT-proBNP standard	Stored at 2-8 °C.
HRP-Conjugate Canine NT-proBNP detection antibody, Standard diluent, Sample Diluent, Chromogen solution A, Chromogen solution B, Stop Solution, Wash buffer	Stored at 2-8 °C.

### PRECAUTIONS

1. Do not substitute reagents from one kit lot to another. Standard, conjugate and microtiter plates are matched for optimal performance. Use only the reagents supplied by manufacturer.
2. Allow kit reagents and materials to reach room temperature (20-25 °C) before

- use. Do not use water baths to thaw samples or reagents.
3. Do not use kit components beyond their expiration date.
  4. Use only deionized or distilled water to dilute reagents.
  5. Do not remove microtiter plate from the storage bag until needed. Unused strips should be stored at 2-8 °C in their pouch with the desiccant provided.
  6. Use fresh disposable pipette tips for each transfer to avoid contamination.
  7. Do not mix acid and sodium hypochlorite solutions.
  8. Serum and plasma should be handled as potentially hazardous and capable of transmitting disease. Disposable gloves must be worn during the assay procedure, since no known test method can offer complete assurance that products derived from Rat blood will not transmit infectious agents. Therefore, all blood derivatives should be considered potentially infectious and good laboratory practices should be followed.
  9. All samples should be disposed of in a manner that will inactivate viruses.
  10. Liquid Waste: Add sodium hypochlorite to a final concentration of 1.0%. The waste should be allowed to stand for a minimum of 30 minutes to inactivate the viruses before disposal.
  11. Chromogen Solution is easily contaminated. If bluish prior to use, do not use.
  12. Chromogen B contain 20% acetone, keep this reagent away from sources of heat or flame.
  13. Remove all kit reagents from refrigerator and allow them to reach room temperature (20-25 °C).

#### **TECHNICAL HINTS**

- When mixing or reconstituting protein solutions, always avoid foaming.
- To avoid cross-contamination, change pipette tips between additions of each

Standard level, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.

- To ensure accurate results, proper adhesion of plate sealers during incubation steps is necessary.
- When using an automated plate washer, adding a 30 second soak period following the addition of Wash Buffer, and/or rotating the plate 180 degrees between wash steps may improve assay precision.
- Chromogen solution should remain colorless until added to the plate. Keep Chromogen solution protected from light. Chromogen solution should change from colorless to gradations of blue.
- Stop Solution should be added to the plate in the same order as the Chromogen solution. The color developed in the wells will turn from blue to yellow upon addition of the Stop Solution. Wells that are green in color indicate that the Stop Solution has not mixed thoroughly with the Chromogen solution.

#### **PRECISION**

Intra-assay Precision (Precision within an assay): Four samples of known concentration were tested twenty times on one plate to assess intra-assay precision. The CV (%) <9%.

Inter-assay Precision (Precision between assays): Three samples of known concentration were tested in twenty separate assays to assess inter-assay precision. Assays were performed by at least three technicians using two lots of components. The CV (%) <11%.

## RECOVERY

The recovery of Canine NT-proBNP spiked to different levels in samples throughout the range of the assay in various matrices was evaluated. The recovery ranged from 98% to 116% with an overall mean recovery of 106%.

## SENSITIVITY

The minimum detectable dose (MDD) of Canine NT-proBNP is typically less than 0.1 pg/ml. The MDD was determined by adding two standard deviations to the mean O.D. value of twenty zero standard replicates and calculating the corresponding concentration.

## SPECIFICITY

Canine N-terminal pro-brain natriuretic peptide ELISA Kit can be used to measure Canine NT-proBNP in samples. Canine N-terminal pro-brain natriuretic peptide ELISA Kit has high sensitivity and excellent specificity for detection of Canine NT-proBNP. No significant cross-reactivity or interference between Canine NT-proBNP and analogues was observed.

## SAMPLE COLLECTION AND STORAGE

The sample collection and storage conditions listed below are intended as general guidelines. Sample stability has not been evaluated.

**Tissue homogenates** - For general information, hemolysis blood may affect the result, so you should rinse the tissues with ice-cold PBS (0.01M, pH=7.4) to remove excess blood thoroughly. Tissue pieces should be weighed and then minced to small pieces which will be homogenized in PBS (the volume depends on the weight of the tissue. 9ml PBS would be appropriate to 1 gram tissue pieces. Some

protease inhibitor is recommended to add into the PBS.) with a glass homogenizer on ice. To further break the cells, you can sonicate the suspension with an ultrasonic cell disrupter or subject it to freeze-thaw cycles. The homogenates are then centrifuged for 5 minutes at 5000×g to get the supernate.

**Cell culture supernatants and other biological fluids** - Centrifuge samples for 20 minutes at 1000×g. Remove particulates and assay immediately or store samples in aliquot at -20 °C or -80 °C for later use. Avoid repeated freeze/thaw cycles.

**Serum** - Use a serum separator tube and allow samples to clot for two hours at room temperature or overnight at 4 °C before centrifugation for 20 minutes at approximately 1000×g. Assay freshly prepared serum immediately or store samples in aliquot at -20 °C or -80 °C for later use. Avoid repeated freeze/thaw cycles.

**Plasma** - Collect plasma using EDTA or heparin as an anticoagulant. Centrifuge samples for 15 minutes at 1000×g at 2-8 °C within 30 minutes of collection. Remove plasma and assay immediately or store samples in aliquot at -20 °C or -80 °C for later use. Avoid repeated freeze/thaw cycles. Samples should be centrifuged adequately and no hemolysis or granule was allowed.

## REAGENT PREPARATION

Bring all reagents to room temperature before use. If crystals have formed in the Buffer Concentrates, warm them gently until they completely dissolved.

**Wash buffer** - Dilute with Distilled or deionized water 1:30.

**Standard** - Pipette 150 µl of Standard Diluent into each tube. Use the stock solution to produce a 2-fold dilution series (below). Mix each tube thoroughly

before the next transfer.


The undiluted Standard serves as the high standard. Standard Diluent serves as the zero standard.

**Note:** If samples generate values higher than the highest standard, please dilute the samples with Sample Diluent and repeat the assay.

### ASSAY PROCEDURE

1. Dilution of standard solutions: This kit contains a standard of known concentration, which could be diluted in small tubes by the end-user by following the instruction in the table below:

60 pg/ml	Standard No.5	150µl Original Standard + 150µl Standard diluents
30 pg/ml	Standard No.4	150µl Standard No.5 + 150µl Standard diluents
15 pg/ml	Standard No.3	150µl Standard No.4 + 150µl Standard diluent
7.5 pg/ml	Standard No.2	150µl Standard No.3 + 150µl Standard diluent
3.75 pg/ml	Standard No.1	150µl Standard No.2 + 150µl Standard diluent

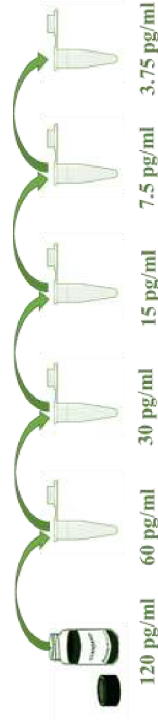

2. Prepare all reagents before starting assay procedure. It is recommended that all Standards and Samples be added in duplicate to the microplate.
3. Add standard: Set Standard wells, testing sample wells. Add diluted standard 50µl to standard well.
4. Add Sample: Add sample diluent 40 µl to testing sample well. Then add 10µl sample to testing sample well; in Blank well don't add anything.

5. Cover with a plate cover and incubate for 45 minutes at 37 °C.
6. Aspirate each well and wash, repeating the process four times for a total of five washes. Wash by filling each well with Wash buffer (250 µl) using a squirt bottle, manifold dispenser or autowasher. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash buffer by aspirating or decanting. Invert the plate and blot it against clean paper towels.

7. Add 50 µl HRP-Conjugate detection antibody to each well, except blank well.

8. Cover with a new adhesive strip. Incubate for 30 minutes at 37 °C.

9. Repeat the aspiration/wash process for five times as in step 5.

10. Add 50 µl chromogen solution A and 50 µl chromogen solution B to each well. Gently mix and incubate for 15 minutes at 37°C. Protect from light.

11. Add 50µl Stop Solution to each well. The color in the wells should change from blue to yellow. If the color in the wells is green or the color change does not appear uniform, gently tap the plate to ensure thorough mixing.

12. Read the Optical Density (O.D.) at 450 nm using a microtiter plate reader within 15 minutes.

## PROTOCOL SUMMARY

Prepare reagents, samples and standards.



Add prepared samples and standards in respective wells together and incubate for 45 minutes at 37 °C.



Wash the plate five times. Add HRP-Conjugate and incubate for 30 minutes at 37 °C.



Wash the plate five times. Add Chromogen solution A and B. Incubate for 15 minutes at 37 °C for color development and add stop solution.



Read the OD value within 15 minutes.



Plot and Calculate.

## CALCULATION

Average the duplicate readings for each standard, control, and sample and subtract the average zero standard optical density (O.D.). Create a standard curve by reducing the data using computer software capable of generating a four parameter logistic (4-PL) curve fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration for each standard on the x-axis and draw a best fit curve through the points on the graph. The data may be linearized by plotting the log of the Canine NT-proBNP concentrations versus the log of the O.D. and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less precise fit of

the data. If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.

**ASSAY RANGE :** 3.75 pg/mL - 60 pg/mL.

**PACKAGE SIZE :** 96Tests.

**SENSITIVITY :** <0.1 pg/ml.

**Linearity:** To assess linearity of the assay, samples containing and/or spiked with high concentrations of Canine NT-proBNP were diluted with the appropriate calibrator diluent to produce samples with values within the dynamic range of the assay. Linear regression analysis of samples versus the expected concentration yielded a correlation coefficient of 0.99.

**VALIDITY & STORAGE:** Twelve months (at 2-8°C, unopened).